

# QuikScreen®

## One Step Onsite Drug Cup

### Available with Adulteration Test

Unlike other onsite tests, QuikScreen® does not require any handling or manipulation to activate the test. Just dispense urine into the cup and read the results. The cup does all the work, IN ONLY 5 MINUTES.



- ★ SIMPLE PROCEDURE
- ★ EASY TO INTERPRET COLOR BANDS
- ★ BUILT-IN TEST CONTROL
- ★ ROOM TEMPERATURE STORAGE
- ★ 18 MONTHS SHELF LIFE
- ★ NO INSTRUMENTS NEEDED
- ★ OPTIONAL ADULTERATION STRIPS BUILT-IN!

QuikScreen® is the easy to use urinalysis drug test that detects drugs and drug metabolites at **SAMHSA** cutoff levels.

#### ORDERING INFORMATION:

	REGULAR CUP	ADULTERATION CUP
6 Test Cup	606XX	656XX
5 Test Cup	605XX	655XX
4 Test Cup	604XX	654XX

"X" denotes number representation of each drug combination. Refer to catalog listing for combination and catalog number.

#### Create your own combination!

Tests available for MET, AMP, MAD, BZD, OPI, THC, COC, BAR, and PCP.

Cups available in five, six, seven, eight, and nine test formats.


# QuikScreen® One Step Onsite Drug Cup

- ★ **Accuracy:** FDA 510(K) registered detection strips have the same high degree of accuracy when compared to SAMHSA certified laboratory screening methods.
- ★ **Self-Contained:** Each cup includes a temperature strip and identification label. Each cup is sealed in its own bag and has an 18-month shelf life.
- ★ **No Leakage:** The cup does not leak so a positive specimen can be shipped in the same container to a laboratory for GC/MS confirmation.
- ★ **Cost:** In addition to providing immediate results, QuikScreen® device costs less than half the price of laboratory testing.



## Test procedure


1. Urinate into the cup, filling the cup with urine above the minimum level line.
2. Wait 5 minutes.
3. Read the results.
4. Refer to product insert for detailed instructions.

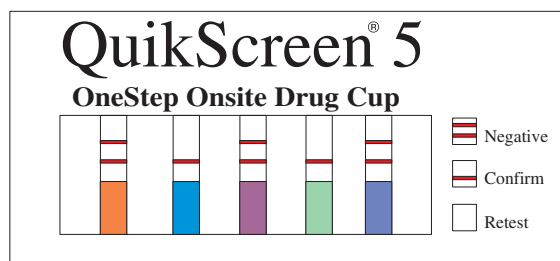

**Confirm:** One line indicates a positive result.

**Negative:** Two lines indicate a negative result.

## Easy to read results

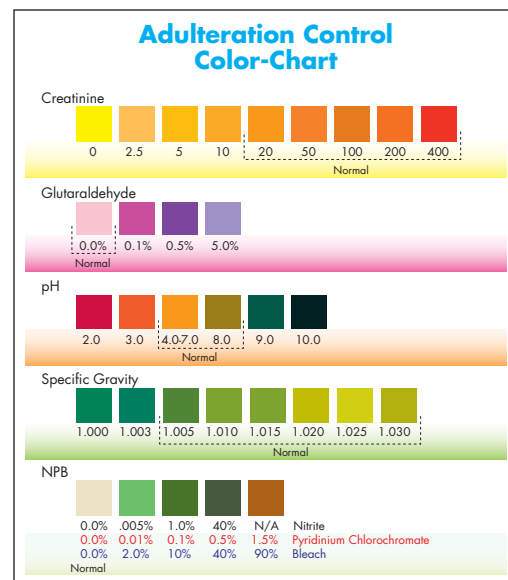

**Compare adulteration results to color chart for devices with built-in adulteration check.**

(Above chart is for illustration purpose only and should not be used to compare results.)

Contact [DrugTesting-kits.com](http://DrugTesting-kits.com) for more information about our QuikScreen® products. Call, fax, or e-mail us your order at [sales@drugtesting-its.com](mailto:sales@drugtesting-its.com). Visit us at [www.drugtesting-its.com](http://www.drugtesting-its.com)

Doc. ID: PDS-057  
Revision: 01